


How to deliver bottles with waste solvents

Never fill in more than 90 % of the total volume!

Collected waste solutions may still react. Then a minimum of 10 % of the volume is needed as an expansion space.

Ensure that there is not reaction potential of your waste bottle:

1. Close it tightly and then agitate it forcefully.
2. Then open the cap a little to enable any developed gas to escape. Leave the bottle standing overnight.
3. Close the cap tightly. Ensure that the bottle is clean and not leaking.

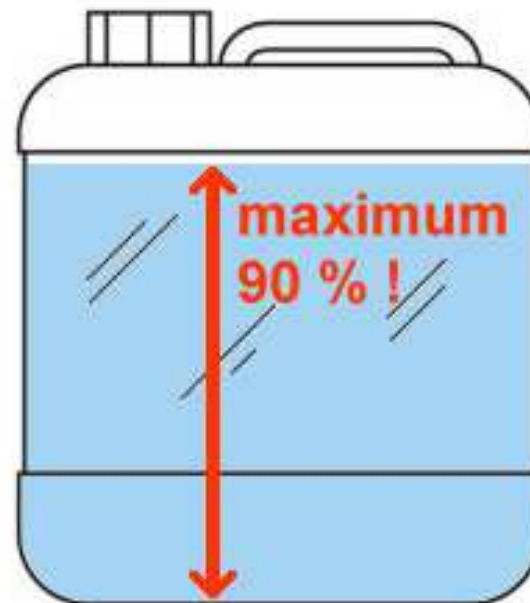

Never put even traces of halogenated compounds to halogen-free waste solvents!